



## Clinical Practice Procedures: Other/Active self-warming blanket – BARRIER® EasyWarm®

<b>Policy code</b>	CPP_OT_ACSB_1221
<b>Date</b>	December, 2021
<b>Purpose</b>	To ensure a consistent procedural approach for the active self-warming blanket – BARRIER® EasyWarm®
<b>Scope</b>	Applies to Queensland Ambulance Service (QAS) clinical staff.
<b>Health care setting</b>	Pre-hospital assessment and treatment.
<b>Population</b>	Applies to all ages unless stated otherwise.
<b>Source of funding</b>	Internal – 100%
<b>Author</b>	Clinical Quality & Patient Safety Unit, QAS
<b>Review date</b>	December, 2024
<b>Information security</b>	UNCLASSIFIED – Queensland Government Information Security Classification Framework.
<b>URL</b>	<a href="https://ambulance.qld.gov.au/clinical.html">https://ambulance.qld.gov.au/clinical.html</a>

While the QAS has attempted to contact all copyright owners, this has not always been possible. The QAS would welcome notification from any copyright holder who has been omitted or incorrectly acknowledged.

All feedback and suggestions are welcome. Please forward to: [Clinical.Guidelines@ambulance.qld.gov.au](mailto:Clinical.Guidelines@ambulance.qld.gov.au)

### Disclaimer

The Digital Clinical Practice Manual is expressly intended for use by appropriately qualified QAS clinicians when performing duties and delivering ambulance services for, and on behalf of, the QAS.

The QAS disclaims, to the maximum extent permitted by law, all responsibility and all liability (including without limitation, liability in negligence) for all expenses, losses, damages and costs incurred for any reason associated with the use of this manual, including the materials within or referred to throughout this document being in any way inaccurate, out of context, incomplete or unavailable.

© State of Queensland (Queensland Ambulance Service) 2021.



This work is licensed under the Creative Commons Attribution-NonCommercial-NoDerivatives V4.0 International License

You are free to copy and communicate the work in its current form for non-commercial purposes, as long as you attribute the State of Queensland, Queensland Ambulance Service and comply with the licence terms. If you alter the work, you may not share or distribute the modified work. To view a copy of this license, visit <http://creativecommons.org/licenses/by-nc-nd/4.0/deed.en>

For copyright permissions beyond the scope of this license please contact: [Clinical.Guidelines@ambulance.qld.gov.au](mailto:Clinical.Guidelines@ambulance.qld.gov.au)

# Active self-warming blanket – BARRIER® EasyWarm®

December, 2021

Effective, continuous patient warming reduces the risks associated with hypothermia. Hypothermia complications include; coagulopathy, delayed wound recovery, decreased oxygen delivery and decreased metabolism.<sup>[1]</sup> Research has shown that hypothermia occurs in 50–60% of all trauma patients.

The **BARRIER® EasyWarm®** is an active self-warming blanket that once opened and unfolded, will heat to 44°C within 30 minutes and maintain a constant temperature for up to 10 hours.



✓ **Indications**

- Hypothermia associated with:
  - trauma
  - environmental exposure

✗ **Contraindications**

- Patients < 18 years
- Use in aeromedical operations

! **Complications**

- The blanket may take up to 30 minutes to warm. Officers should ensure the blanket is opened and prepared as soon as possible.
- The blanket must not be positioned under the patient.
- Direct blanket to skin contact in areas with bruising, swelling and frostbite should be avoided.
- In patients with impaired sensibility, reactivity or communication, skin response should be regularly monitored for signs of over exposure.
- The blanket must not be folded over itself.
- Officers should avoid covering pads with belts or fixation straps.

## Procedure – Active self-warming blanket – BARRIER® EasyWarm® [2]

1. Tear package open using pre-cut located on the top left corner of the packaging.
2. Remove blanket from the protective package and expose to air.
3. Open and completely unfold the blanket.


Pre-cut tab  
(tear to open)



4. Place the blanket on the patient with the pads facing up.
5. Consider placing one (1) open weave cotton blanket over the BARRIER® EasyWarm® for additional coverage.
6. Regularly monitor the patient for signs of redness indicating over exposure.

Air activated **warming pads** integrated into pocket on the top layer.

The warming pads are made up of air permeable bags containing active coal, salt, water and iron powder.

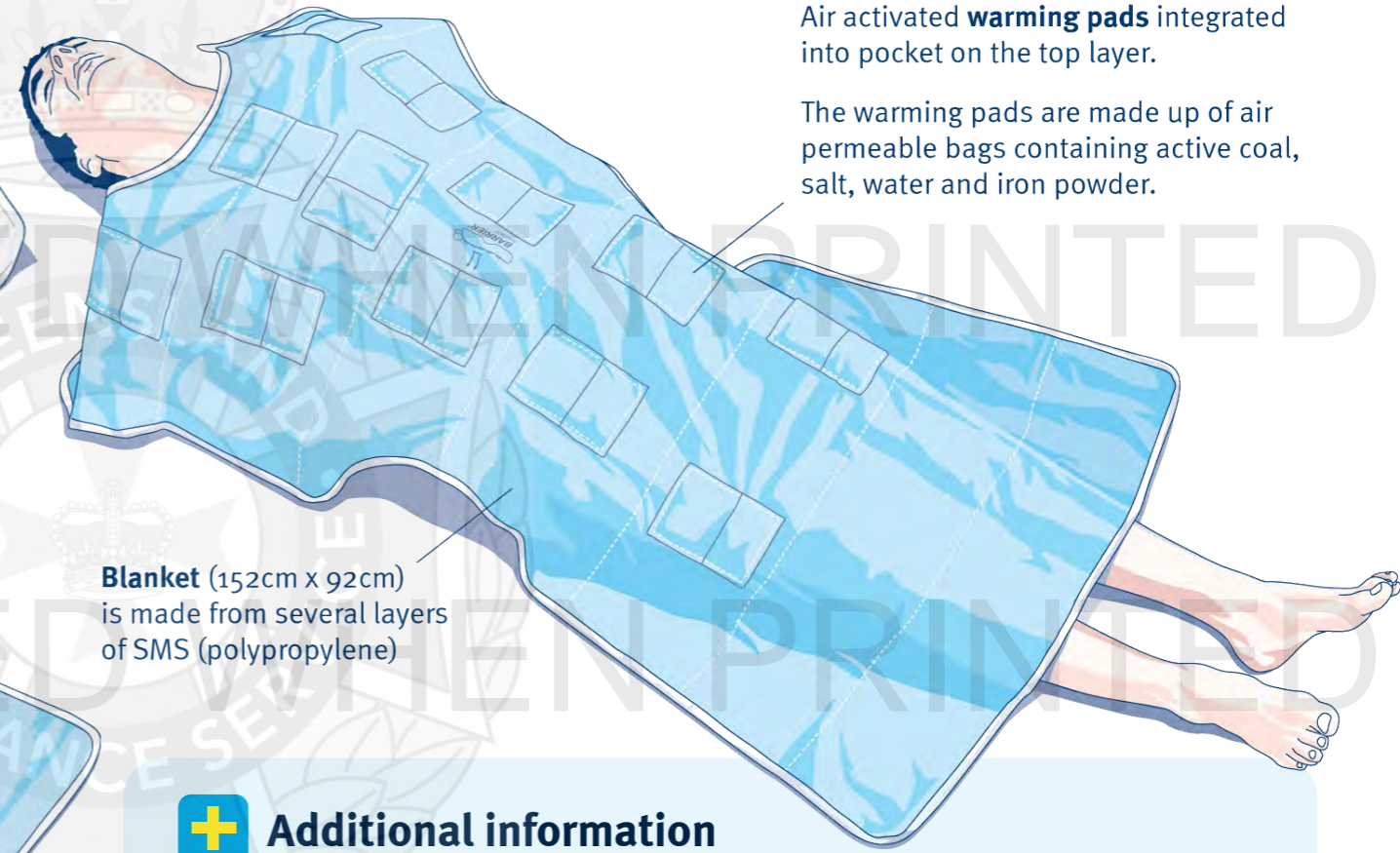

**Blanket** (152cm x 92cm)  
is made from several layers  
of SMS (polypropylene)

### + Additional information

- Officers should consider other warming measures including removal of wet clothes and in-vehicle heating prior to application of the BARRIER® EasyWarm®.
- On arrival at the health facility clinicians must ensure that responsibility for ongoing monitoring of the patient's skin response is handed over to the hospital staff.

